Achtung:

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Attention:

Please do not dispose of old electronic devices and empty batteries in household waste. To protect the environment, take them to your retail store or to appropriate collection sites according to national or local regulations.

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- ① ON/OFF button
- ② Sensor
- ③ Display window

2. FEATURES

- LCD displays the measured temperature.
- With memory for the last 5 measurements
- With beeper to alarm the peak temperature.
- Auto self-test and power auto shut-off.
- Safe – without glass and mercury.

3. SPECIFICATIONS

Temperature Range:	32.0 ... 43.0 °C
Accuracy:	±0.1°C 35.0 ... 41.0 °C ±0.2°C the rest
Resolution:	0.1 °C
Working Conditions:	0 ... 50°C Humidity: 10 ... 83%
Storage Conditions:	-20 ... 60°C Humidity: 10 ... 83%
Battery:	1.55Vdc Alkaline Manganese battery SR41, LR41, UCC392, UCC 192
Battery Life:	approx. 100 hours for continuous operation or approx. 18 months while used 10 minutes per day
Power consumption:	0.15mW
Sensor Type:	Thermistor
Display:	Liquid Crystal Display 3-1/2 digits
Memory:	Memory function for the last 5 measured temperature values
Alarm:	approx. 8 seconds of alarm sound while peak temperature is stable
Dimension:	55 x 50 x 50 mm
Weight:	approx. 12 grams with battery inside

4. HOW TO USE

- Press ON/OFF button to switch on. An acoustic signal occurs and the LCD displays 188.8, and then 36.5°C followed by L°C (if the ambient temperature is lower than 32.0°C). When the °C symbol starts flashing, the thermometer is ready for use.
- To recall the last 5 readings, press and hold ON/OFF button for 2 sec. while switching on. The LCD displays 188.8, followed by 1M and the last measured value appears. Press ON/OFF button again to recall the preceding values 2M-5M.
- If the ambient temperature is higher than 43.0°C, the LCD displays H°C, if it is lower than 32.0°C it displays L°C.
- Position the pacifier bulb in child's mouth and allow the child to suck on the Pacifier Thermometer for a minimum of 3 minutes (allowing the thermometer to remain in the child's mouth for 6 to 8 minutes will provide the most accurate reading).
- If the peak temperature is stable for 32 seconds, the alarm sound occurs for 8 seconds.
- Remove the Pacifier Thermometer and read the temperature indicated on LCD.
- To turn the thermometer off, press the On/Off button again. In order to conserve battery life, we recommend turning the thermometer off after taking a reading. However, the thermometer will automatically shut off after 10 minutes.
- If the measurement is interrupted for any reason, repeat the above procedure.

5. PRECAUTIONS

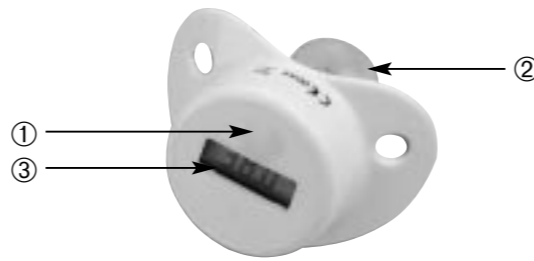
- Examine the pacifier bulb before each use for cracks, tears, tackiness or other signs of wear. Discard the worn unit immediately.
- Do not expose to direct sunlight, high temperature moisture or dirt.
- Do not attempt to pull or twist the nipple apart from the body or detach button bead as this will damage the Thermometer's sensitive temperature probe and electronics.
- Never allow a child to play with this product. It is not a toy. It is a sensitive electronic device and maybe damaged if not used properly under parental supervision.
- Use this Thermometer for oral temperature measurements only.
- The battery is not replaceable. Do not attempt to remove the battery.

6. CLEANING

- The thermometer may be sterilized by immersing the silicone rubber nipple ONLY in sterilizing solution such as alcohol. Immerse the nipple approximately 1 second ONLY in the solution. The protective clip-on plastic shield supplied maybe used as container for this purpose.
- **Warning:** The Pacifier Thermometer should never be completely immersed in sterilizing solution or sterilized in boiling water as damage will occur to the electronic components.

7. LIMITED WARRANTY

- Contact the original point of purchase for claiming. No guarantee if the instrument is handled or opened improperly.
TFA Dostmann GmbH & Co. KG · 97877 Wertheim-Reicholzheim · Germany

**(GB)****Bedienungsanleitung
Operating Instructions****Mode d'emploi****Istruzioni****Gebruiksaanwijzing****Instrucciones para el manejo****CE****Kat. Nr. 15.2005****Thermomètre-tétine**

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